



# Pevco Tracker™

## Real-Time Delivery Confirmation

- origin
- time sent
- contents
- destination
- time arrived



# Pevco Tracker™

Pevco Tracker™ puts **real-time pneumatic tube system transaction data** and arrival confirmations directly into the hands of the **healthcare professional** responsible for providing timely patient care.



Pevco Tracker™ is ideal for Pharmacies and Blood Banks and is available with an optional antimicrobial housing.

## Improve Patient Care and Reduce Errors

### Confirm Deliveries

Lab and pharmacy personnel receive status messages via a Tracker PDA when it is important to know exactly when a pharmaceutical, blood product or lab specimen has arrived at its destination.


### Avoid Delays

Misplaced items can be traced and located before they perish by reviewing detailed transaction logs on a PDA or networked computer.

### Improve Performance

Managers can analyze detailed historical transaction records stored in Tracker's database to gain insight into department activity and performance.

Automated messages provide delivery confirmation, as well as alert staff of incoming carriers or possible delays.

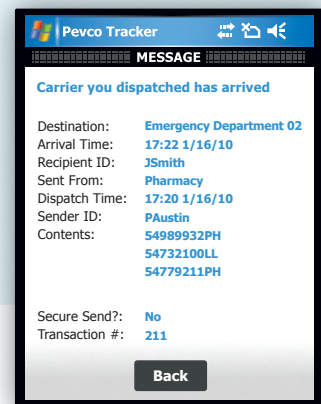

## > Easy to Implement and Expand

Installed with little or no system down-time, hospitals can use Pevco Tracker™ at every pneumatic tube station or only at a few where needed. PDA's can be added to the system as staffing or transaction needs change. Pevco Tracker™ relies on proven barcode technology and merges seamlessly with your existing Pevco pneumatic tube system. Equipped with familiar touch screen operations and clear, easy to read graphics, Tracker is simple for staff to learn and use.

## > Flexible

Pevco Tracker™ provides administrators with the flexibility to choose which elements of a transaction need to be captured. Senders may choose to scan only the contents and the destination address. But when higher accountability is needed, every element can be scanned recording the "who, what, when and where" of each transaction. For deliveries which require higher security, a secure transaction option is available directly from the touch screen. In a secure send, the incoming carrier is held above the receiving station until the recipient scans their badge ensuring delivery documentation and safe arrival.

## > Reduces Errors

With Pevco Tracker™, healthcare professionals have an easy to use tool that increases visibility and accountability of all pneumatic tube system transactions. High-value medications, specimens and blood products can be tracked from sender to receiver. Carriers sent to the wrong destination because of user error can be tracked instantly from a Tracker PDA or a networked PC and quickly rerouted to reduce product loss and ensure timely patient care. Optional PDA alert messages offer staff time to react to errors quickly.



## Pevco Tracker™ supports both standard and secure pneumatic tube system transactions:

### Standard Transaction

- What was sent
- When it was sent
- When it arrived
- Where it was sent from
- Where it was sent to

### Secure Transactions

- What was sent
- Who sent it
- When it was sent
- When it arrived
- Who received it
- Where it was sent from
- Where it was sent to

For more information on Pevco Tracker™ and other reliable Pevco products visit [www.pevco.com](http://www.pevco.com) or call 1.800.296.7382.



Healthcare's Pneumatic Tube Delivery System Expert™

[www.pevco.com](http://www.pevco.com)

**Pevco™**

1401 Tangier Drive, Baltimore, Maryland 21220  
Tel: 410.931.8800 800.296.7382 Fax: 410.931.4660